# RadiForce® MX270W

## Color LCD Monitor

### Instructions for Use: ■ Setup Manual

**Important:** Please read PRECAUTIONS, this Setup Manual and the User's Manual stored on the CD-ROM carefully to familiarize yourself with safe and effective usage. Please retain this manual for future reference.

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## Compatible Resolutions

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The monitor supports the following resolutions.

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
# RadiForce® MX270W

Color LCD Monitor

## Location of the Caution Statements
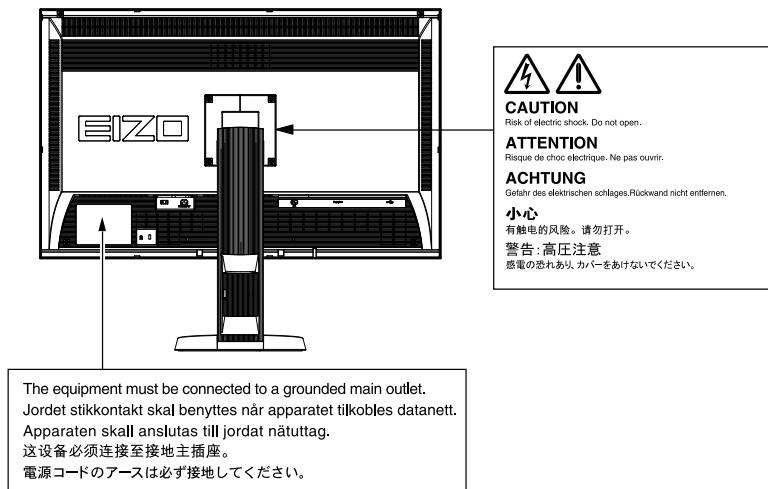
## Symbols on the unit

Symbol	This symbol indicates	
	Main Power Switch:	Press to turn the monitor's main power off.
	Main Power Switch:	Press to turn the monitor's main power on.
	Power button:	Press to turn the monitor's power on or off.
	Alternating current	
	Alerting electrical hazard	
	CAUTION:	Refer to SAFETY SYMBOLS section in PRECAUTIONS.

No part of this manual may be reproduced, stored in a retrieval system, or transmitted, in any form or by any means, electronic, mechanical, or otherwise, without the prior written permission of EIZO Corporation. EIZO Corporation is under no obligation to hold any submitted material or information confidential unless prior arrangements are made pursuant to EIZO Corporation's receipt of said information. Although every effort has been made to ensure that this manual provides up-to-date information, please note that EIZO monitor specifications are subject to change without notice.


## Package Contents

Check that all the following items are included in the packaging box. If any items are missing or damaged, contact your local dealer.


### Note

- Please keep the packaging box and materials for future movement or transport of the monitor.


- Monitor
- Power cord

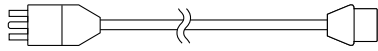
- Digital signal cable (DisplayPort) : PP300

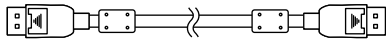- Digital signal cable (Dual Link) : DD300DL

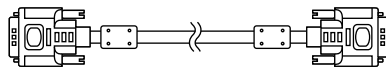
- EIZO USB cable: UU300

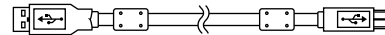- EIZO LCD Utility Disk (CD-ROM)
- Setup Manual (this manual)
- PRECAUTIONS

## Controls and Functions
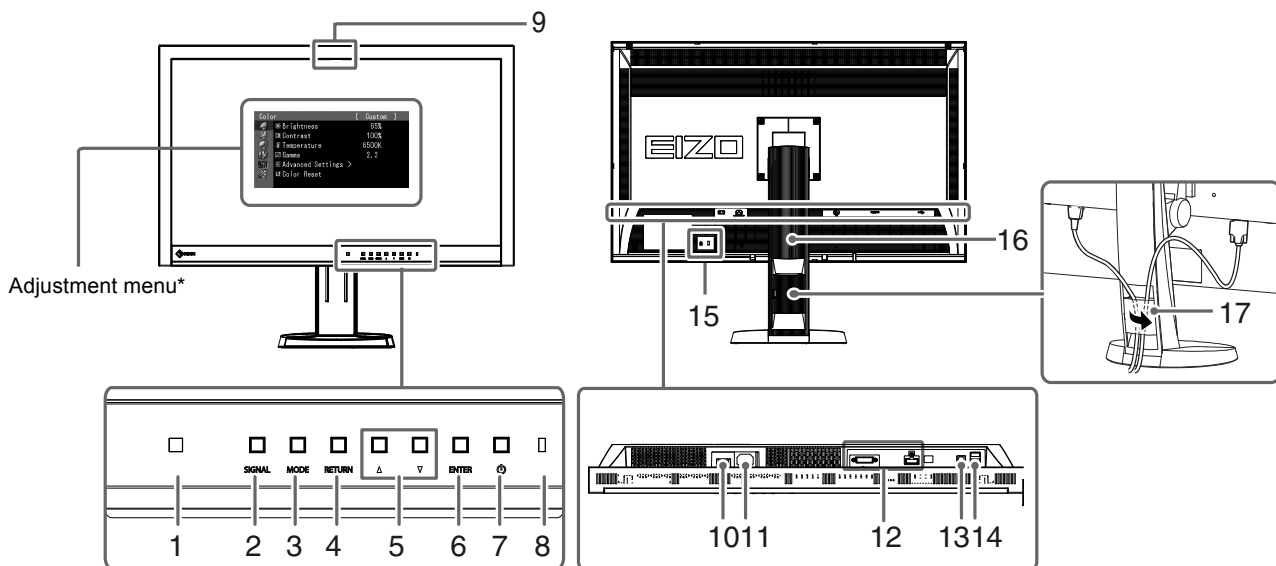
<b>1. Ambient light sensor</b>	Used for the illuminance change monitoring function of RadiCS/RadiCS LE. For details, refer to RadiCS/RadiCS LE User's Manual (on the CD-ROM).
<b>2. SIGNAL button</b>	Switches input signals for display.
<b>3. MODE button</b>	Switches the CAL Switch mode.
<b>4. RETURN button</b>	Cancels the setting/adjustment and exits the Adjustment menu.
<b>5. ▲▼ button</b>	Chooses and increases/decreases adjustment settings for advanced adjustment using the Adjustment menu.
<b>6. ENTER button</b>	Displays the Adjustment menu, determines an item on the menu screen, and saves values adjusted.
<b>7. ⏻ button</b>	Turns the power on or off.
<b>8. Power indicator</b>	Indicates the monitor's operation status. Blue: Operating, Orange: Power saving mode, Off: Main power/Power off
<b>9. Integrated Front Sensor (slide type)</b>	Used for calibration and Grayscale Check.
<b>10. Main power switch</b>	Turns the main power on or off.
<b>11. Power connector</b>	Connects the power cord.
<b>12. Input signal connector</b>	Left: DVI-D connector / Right: DisplayPort connector
<b>13. USB upstream port</b>	Connects the USB cable to use the software that needs USB connection, or to use USB Hub function.
<b>14. USB downstream port</b>	Connects a peripheral USB device.
<b>15. Security lock slot</b>	Complies with Kensington's MicroSaver security system.
<b>16. Stand</b>	Used to adjust the height and angle of the monitor screen.
<b>17. Cable holder</b>	Covers the monitor cables.

\* For instructions for use, refer to the User's Manual (on the CD-ROM).

# Connecting Cables

## Attention

- Check that the monitor and the PC are powered off.
- When replacing the current monitor with this monitor, be sure to change the PC settings for resolution and vertical scan frequency to those that are available for this monitor. Refer to the compatible resolution table (back of cover page) before connecting the PC.


## Note

- When connecting multiple PCs to this product, refer to the User's Manual (on the CD-ROM).

## 1. Connect the signal cables to the input signal connectors and PC.


Check the shape of the connectors, and connect the cables.

After connecting the signal cable, tighten the screws of the connectors to secure the coupling.
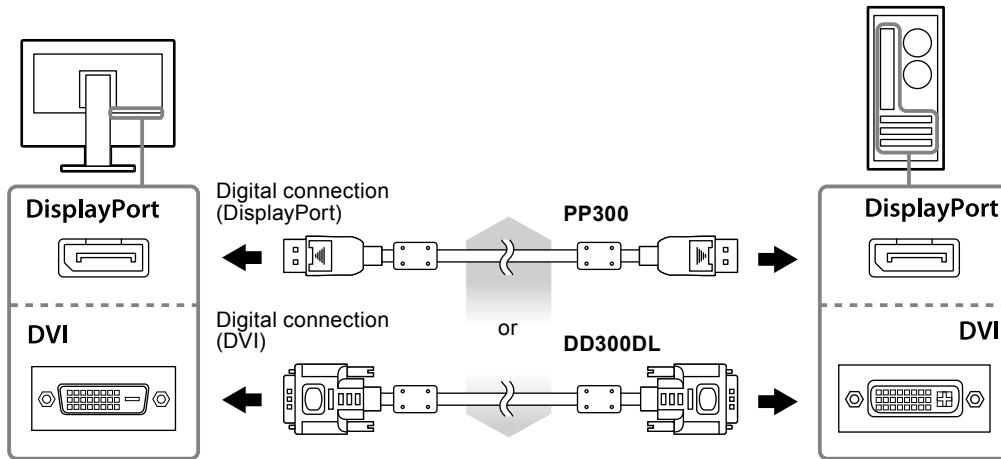
## 2. Plug the power cord into a power outlet and the power connector on the monitor.

## 3. Connect the USB cable when using RadiCS LE or ScreenManager Pro for Medical.
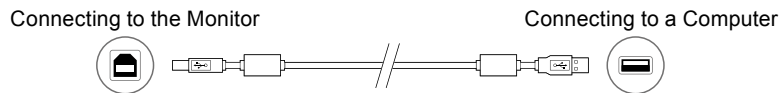
## 4. Press to turn on the monitor.

The monitor's power indicator lights up blue.

## 5. Turn on the PC.

The screen image appears.


If an image does not appear, refer to "No picture" (page 6) for additional advice.

## Attention

- Turn off the monitor and PC after using them.
- For the maximum power saving, it is recommended that the Power button be turned off. Turning off the main power switch or unplugging the power cord completely shuts off power supply to the monitor.

# Adjusting the Screen Height and Angle


Hold left and right edge of the monitor with both hands, and adjust the screen height, tilt and swivel of the screen to the best condition for working.

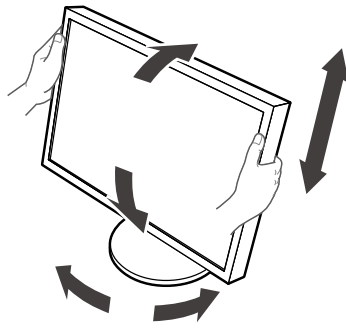
## Attention

- Be sure that the cables are correctly connected.

# Showing Button Guide

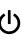Press the front buttons (except  button), the Button Guide appears above the button.

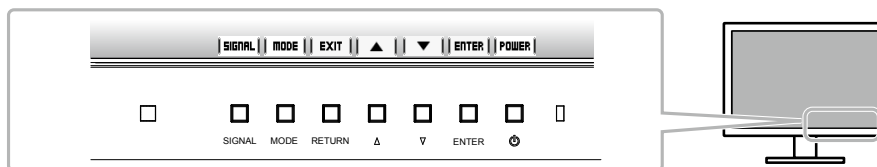
## Note

- The Button Guide will continue to appear while the Adjustment menu or mode menu is showing.
- The Button Guide is displayed differently depending on the selected menu or status.

# Selecting Display Mode

This function allows easy selection of a desired mode according to monitor application.

Mode	Purpose
1-DICOM	Available to display in DICOM mode.
2-Custom	Available for the setting according to your preference.
3-CAL	Displays the screen adjusted by calibration software.
4-Text	Suitable for displaying texts for word processing or spreadsheets.
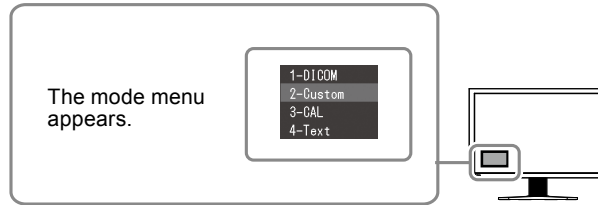
## Procedure

1. Press MODE.  
The mode menu appears at the bottom left of the screen.
2. The mode among the list is highlighted in turn each time MODE is pressed.  
You can switch the mode with ▲ or ▼ while the mode menu is displayed.

# Controlling Monitor Quality

The quality control software “RadiCS LE” (for Windows) is included with this monitor.

To calibrate the monitor and manage the calibration history, install RadiCS LE with reference to the RadiCS LE User's Manual (on the CD-ROM). Use the RadiCS SelfQC function and built-in Integrated Front Sensor to perform periodic calibration and Grayscale Check for the monitor independently. For details, refer to the User's Manual (on the CD-ROM).

### Attention

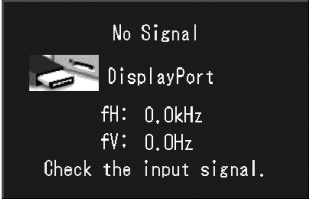
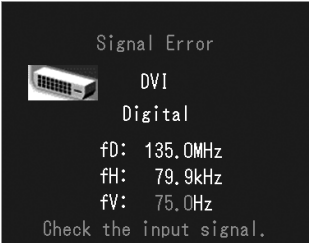
- Connect the monitor to the PC with a USB cable when using RadiCS LE.

### Note

- To perform advanced quality control that conforms to medical standards, use the optional monitor quality control software kit “RadiCS UX1”.

# No picture

If a problem still remains after applying the suggested remedies, contact your local dealer.

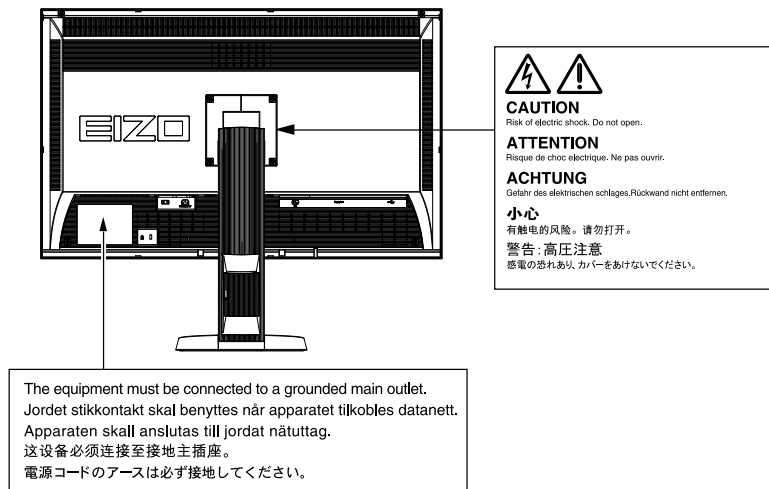
Problem	Possible cause and remedy
<p><b>1. No picture</b></p> <ul style="list-style-type: none"> <li>Power indicator does not light.</li> </ul> <hr/> <ul style="list-style-type: none"> <li>Power indicator is lighting blue.</li> </ul> <hr/> <ul style="list-style-type: none"> <li>Power indicator is lighting orange.</li> </ul> <hr/> <ul style="list-style-type: none"> <li>Power indicator flashes orange and blue.</li> </ul>	<ul style="list-style-type: none"> <li>Check whether the power cord is connected properly.</li> <li>Turn the main power switch on.</li> <li>Press .</li> <li>Turn off the main power, and then turn it on again a few minutes later.</li> </ul> <hr/> <ul style="list-style-type: none"> <li>Increase "Brightness", "Contrast", and/or "Gain" in the Adjustment menu.</li> </ul> <hr/> <ul style="list-style-type: none"> <li>Switch the input signal with SIGNAL .</li> <li>Operate the mouse or keyboard.</li> <li>Check whether the PC is turned on.</li> </ul> <hr/> <ul style="list-style-type: none"> <li>There is a problem in the device that uses DisplayPort connection. Solve the problem, turn off the monitor, and then turn it on again. Refer to the user's manual of the device that uses DisplayPort connection for details.</li> </ul>
<p><b>2. The message below appears.</b></p> <ul style="list-style-type: none"> <li>This message appears when no signal is input. Example:</li> </ul> <div data-bbox="317 887 627 1084" data-label="Image">  </div> <ul style="list-style-type: none"> <li>The message shows that the input signal is out of the specified frequency range. Example:</li> </ul> <div data-bbox="317 1198 627 1440" data-label="Image">  </div> <p>fD: Dot Clock fH: Horizontal Scan Frequency fV: Vertical scan frequency</p>	<p>This message appears when the signal is not input correctly even when the monitor functions properly.</p> <ul style="list-style-type: none"> <li>The message shown left may appear, because some PCs do not output the signal soon after power-on.</li> <li>Check whether the PC is turned on.</li> <li>Check whether the signal cable is connected properly.</li> <li>Switch the input signal with SIGNAL .</li> </ul> <hr/> <ul style="list-style-type: none"> <li>Check whether the PC is configured to meet the resolution and vertical scan frequency requirements of the monitor (see "Compatible Resolutions").</li> <li>Reboot the PC.</li> <li>Select an appropriate display mode using the graphics board's utility. Refer to the manual of the graphics board for details.</li> </ul>

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